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1. LCD DISPLAY
2. MEASUREMENT AREA OF BODY FAT
3. BASE CRYSTAL 2.2CM.
4. LOWER SUPPORT POINT
5. STICKED LEGS
6. BATTERY COMPARTMENT

Operating Instructions

1. SELF-CALIBRATION (TRACK THE ZERO POINT)

The scale is step on technology, performing calibration is a must while using the scale for the first time, or when moving its position or location, they should be restarted. Activate the scale and waiting about 10 seconds until it automatic deactivates on its own (means waiting until the "0,0" disappear from the LCD screen). Then you can stand on the platform directly for weighing. This is for calibrating the scale on a firm, flat surface to let the scale return to ZERO. If you do not perform these measures, the first weight may be incorrect and should be disregarded.


2. WARNING

- People with pacemakers or other types of electronic implants should not use this product because the measuring method of the resistance might interfere with its good function.
- Pregnant woman is not suggested to use this product for body analysis.
- This product is suitable for adults and children of age 10 or older.


- Ensure scale on a flat and firm surface, Avoid carpet or soft surfaces.

3. POWER SUPPLY

This appliance works with battery of 2 pieces AAA ,1.5V. For the first time use, please remove the insulation sheet from the battery jar. If need to change new battery, Install the battery under correct positive and negative plate.

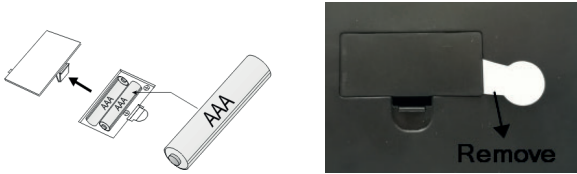4. SCALE OUTLINE

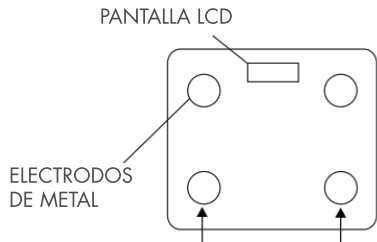Make sure to take off your shoes and socks and make good contact between bare feet and the metal electrodes

5. TECHNICAL PARAMETERS

- 5-1. Turn on: Step On the platform to power on
- 5-2. Auto power off : the scale will turn off automatically after 10 seconds of no operation
- 5-3. Maximum Capacity: 180 kg
- 5-4. Minimum Weight: 5kg
- 5-5. Graduation: 0.1kg
- 5-6. Overload indication: Err
- 5-7. Low battery indication: Lo
- 5-8. Impedance error :Err1 , review and follow the measures in point 14
- 5-9. Memory Person Number: Automatic identify 8 persons
- 5-10. (fat) Measuring body fat range: 5%~50%
- 5-11. (TBW) Measuring water range: 10%~80%

CHECK YOUR MOBILE DEVICE MODEL

Bluetooth 4.0+ and IOS7.0+

Phone 4s	iPhone 5	iPhone 5c
iPhone 5s	iPhone 6	iPhone 6 Plus
iPod touch 5	iPad 3	iPad 4
iPad Air	iPad Air 2	iPad mini
iPad mini 2	iPad mini 3	Or more style

Bluetooth 4.0+ and Android4.3+

Samsung: Galaxy s7	Galaxy s7 edge	Galaxy Note5 etc.
LG LG G5	LG G4	LG V10 etc.
Sony: Xperia Z5	Xperia Z3	Xperia Z2 etc.

And other support for Bluetooth 4.0+ and Android 4.3+ phone or tablet PC models.

APP INSTALATION

1. Turn on Bluetooth
2. The phone or tablet is connected to the network (WIFI, 4G or 3G).
3. Use your device to scan the QR code or "OKOK International" from the APP Store or Google Play to install the APP.


If Android or IOS cannot download the "OKOK International" directly through scanning the barcode, you can go to "Google play" to download it.


USE THE DEVICES AND REGISTER, FILL IN THE USER INFORMATION

Open "OKOK International".


1. If you have the account, fill in the email and password to log in.
2. If you do not have the account, click "Join us", register your personal information correctly. Nickname>Gender>Care>Height>Objective(Target: According to your physical information, your standard weight range is calculated using BMI)
3. If you have registered and have forgotten your password, select "Forgotten password".


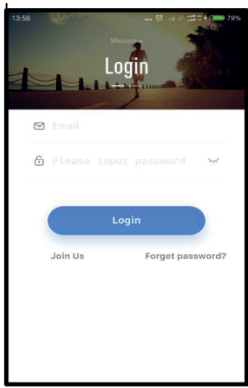
Configure the phone open "bluetooth"; Enter the main APP screen

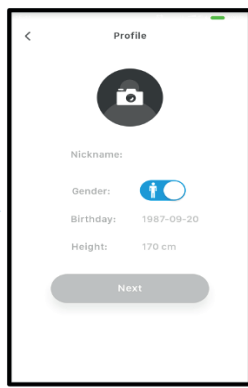Click on "Device"; Switch on the scale; Connect the bluetooth scale to the APP

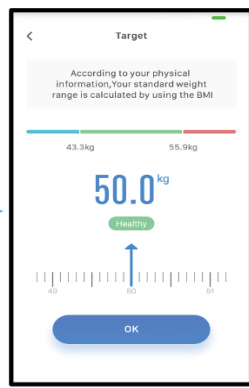Click the "Go" button and add the device.


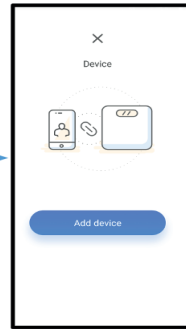
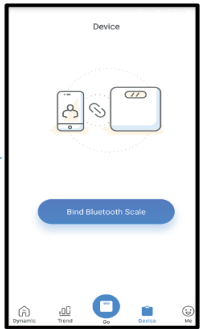
MEASUREMENT

Keep the bluetooth on and stand on the scale, the value will be shown on the scale display.

Click on the recorded value to view the measurement results.

Click on the "Trend" to check the data for each item.

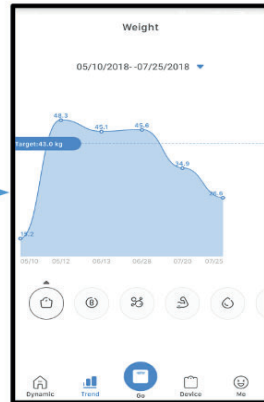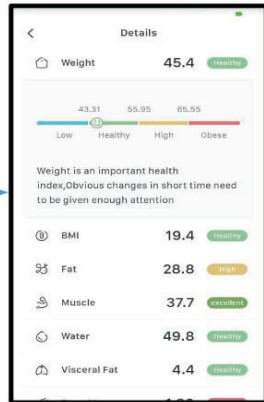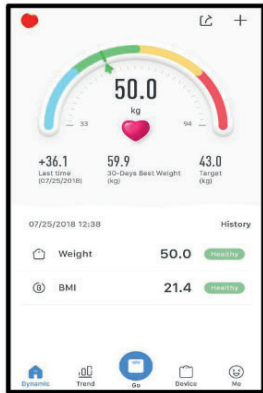APP FUNCTIONS

Click on the picture icon to go to the personal profile, you can change/monitor personal information accordingly.


Displays weight, body fat or BMI, you can scroll down to view histories. Click on the "weight" area to view your weight, fat (%), muscle (%), hydration (%), BMI, visceral fat (%), bone (LB), calories (Kcal/day). Your personal information is related to the accuracy of the algorithm, so make sure your personal data, such as gender, height and age, are entered correctly.

Note:

- 1) If only weight data is available, clicking "weight" will simply display the weight and BMI data..
- 2) The report is based on the complete evaluation of all measurements, so it will only work when all measurements are valid, and all data is automatically saved while weighing.
- 3) If your phone is not available for pairing when you are weighing, the data will be updated when connected.

It will show you the maximum, minimum and average weight, body fat percentage, hydration, muscle mass percentage, bone mass percentage in weekly, monthly and quarterly periods. generated graphs help you to easily see the trend of this data.

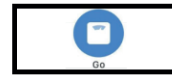Open the APP "OKOK", click "GO" to add the tip before weight, otherwise you can not display the same data as the scale.

1. Open the bluetooth function on the phone, turn on the scale, keep the display light on, link the bluetooth scale.

2. If a scale has connected "OKOK" APP by phone, now you want to change another APP of phone use to connect the scale, firstly, Unlink the bluetooth scale.


1.  My target weight Edit the target weight here.
2.  Family Member Invite 7 other family members and add their name, gender, date of birth and height.
3. The setting interface is used to set the count and unlink the scale and so on.

PROBLEM	CAUSE	SOLUTION
The application does not connect to the scale.	There is a Bluetooth Connection problem between the scale and mobile phone.	<ul style="list-style-type: none"> •Make sure your Bluetooth is active. •Make sure your scale has a battery installed and works. •Switch all other Bluetooth devices nearby off. •Remember that one device can be connected to the scale only. •Close the application and Bluetooth in other devices which use the scale
The scale doesn't activate	<ul style="list-style-type: none"> • No battery. • Battery inserted wrongly • Low battery • Weight below 5 kg. 	<ul style="list-style-type: none"> • Install battery. • Remove battery and insert again, pay attention to battery terminals. • Replace battery. • Minimum measured value is 5 kg.
Serious measurement Error	<ul style="list-style-type: none"> •Scale is not positioned on flat surface. • Check selected measurement unit. • When measuring for the first time, the scale could have not calibrated automatically 	<ul style="list-style-type: none"> • Scale must be positioned on a flat and hard surface. Avoid carpets and other soft surfaces. • Set a proper unit. <p>Please follow up POINT 1 in page 1 for the self-calibration measure.</p>
The display shows Lo	Low battery.	Replace battery.
The display Shows Err	Maximum permissible value has been exceeded.	Maximum permissible value is 200 kg.
The display shows Err1	Scale has difficulties In measuring through electrodes	<ul style="list-style-type: none"> •Measure yourself when you are barefoot, remove footwear and socks •Moisten your feet slightly if necessary. •Remove any thickened areas of the epidermis on your foot. •Make sure you have entered correct values in the application profile (gender, age, height).

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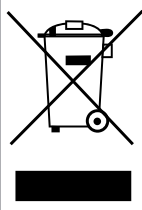INFORMATION FOR USERS.

According to the European guidelines 2002 / 95e, Of the use of dangerous substances in electrical appliances, in addition to waste disposal. The crossed-out symbol of the container in the box indicates that the product at the end of its useful life must deposit in a place separated from the rest. Therefore, the user of delivering the device, when it is no longer used, to the selective collection center for electronic waste, You must return it to the seller at the time of purchase of a new attachment of the same type, one or another move. The separate collection of the unused unit for further recycling, treatment and waste compatible with the environment helps to avoid possible negative effects.